


The Infusion Company

Technical Data Sheet

Product Name	Syramed Extension Line		
Product Code	6210003021 (APEP100M)		
Product Description	Syramed Extension, Polyethylene, ø 1 x 2, 100cm		
Features	 PE line  100 cm  1 ml (40°C) / 1 ml (40°C - 2 bar)  internal Ø: 1.0 mm / outer Ø: 2.0 mm	 Set for pressure equipment  PVC free line	
Product Class	Disposable medical device – IIa Class		
Intended Use	Extension line for syringe used for physiological (saline) infusion and pharmaceutical solutions		
CND Code	A03020101		
Single Packaging	Blister (medical paper + transparent plastic film)	UDI-DI	08052747954506
Multiple Packaging	(24x39x39) cm Box - Qty: 200.0 pcs	UDI-DI	18052747954503
Product Duration	5 years		
Storage Method			
Method of Sterilization	Ethylene Oxide – Residual ETO ≤ 10ppm (ref. “Farmacopea Europea Ufficiale”)		
Disposal Method	The device must be disposed according to local policy for hospital waste management		
Standards	Compliant to ISO 8536-9, ISO 80369-7, ISO 10993-1, ISO 11135, ISO 11737-1, ISO 11607-1, ISO 11607-2		
EC Certificate	0425-MDR-030019-00		
Manufacturer	Phoenix s.r.l Via Leonardo da Vinci - 41038 San Felice sul Panaro (MO) mail: info@phoenixbiomed.it		
Distributor	arcomed ag Steinackerstrasse 29 8302 Kloten, Switzerland www.arcomed.com		

Technical Drawing

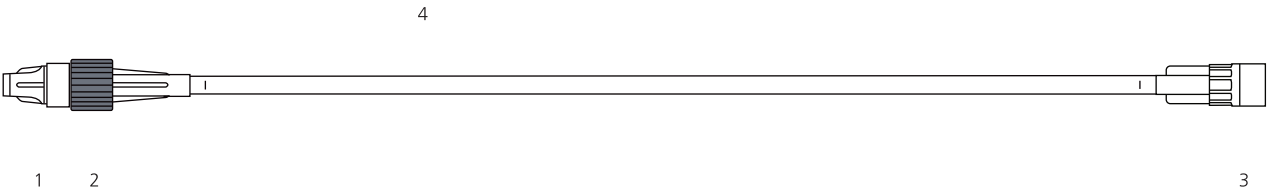

Nr.	Components list: 6210003021
-----	-----------------------------

1	Protection Cap
2	Turnable Ring (blue)
3	Protection Cap
4	ø 1 x 2 Tubing Polyethylene 100 cm - Female-Male Luer Lock Connectors

Warnings

Indications & Notes

Do not use with light-sensitive infusion solutions



Do not re-use

Do not use with Injection Pressure > 30 psig



Do not use if package is damaged



Do not re-sterilize



Caution is necessary when operating the device